



EXray[®] 30

PEEK headless
cannulated compressive
screws
for small fragments



Ø 3,0 mm



SynchroMedical

The interphalangeal reference

EXray® 30 PEEK

Self compressive cannulated screw

DESCRIPTION

EXRAY® 30 compression screws have a core diameter of 2 mm and an outer diameter of 3 mm. The screw is prepared with a 2 mm drill bit and subsequently with a tap to optimize its proper anchorage.

The EXRAY® 30 compression screw is available in sizes from 13 mm to 31 mm.

EXRAY® 30 screws are manufactured in PEEK as per standard ASTM F2026.


- Ø3,0mm Peek cannulated headless compressive screws
- PEEK screws are radiolucent
- The elasticity modulus is close to bone
- Drilling and tapping are mandatory
- Preparation first cortical is mandatory
- Hexalobe T7



ORDERING INFORMATION

Implants

HVP3011.....L 11 mm HVP3023.....L 23 mm
 HVP3013.....L 13 mm HVP3025.....L 25 mm
 HVP3015.....L 15 mm HVP3027.....L 27 mm
 HVP3017.....L 17 mm HVP3029.....L 29 mm
 HVP3019.....L 19 mm HVP3031.....L 31 mm
 HVP3021.....L 21 mm

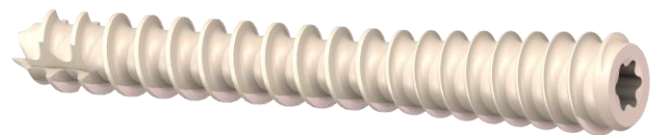

Instruments

	HV9010	AO HANDLE
	HV9140	LONG CANNULATED T7 SCREWDRIVER
	HV9005	GAUGE LENGTH
	HV9011-1	K-WIRE CASE
	HV9016	SCREW EXTRACTOR
	HV9131-013	CANNULATED TAP
	HV9131-031	
	KW10100TR	KIRSCHNER WIRE WITH TROCAR TIP
	HV9040-013	CANNULATED STEP DRILL
	HV9040-031	
	HV9115	EXRAY® INSTRUMENT TRAY
	HV9020	DEPTH GAUGE

All rights reserved © Pictures non-contractual